
# 

## MGUARD BARRIER LEVEL 2 FACEMASK

Enclosed aluminium nose closure forms a strong seal around nose and cheekbones for maximum protection.

The elastic, non-latex ear loops minimise pressure to the ears and are \_\_\_\_\_ attached to the outside of the mask to provide a superior facial seal & fit.

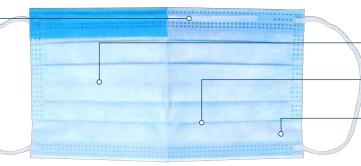

Fluid resistant, spunbond outer layer for maximum protection.


Highest quality meltblown filter media for best filtration efficiency and breathability.

Sonic bonding ensures no pin holes & layer separation.

### OMGuard Barrier Level 2 Mask Earloop 120

#### Product Code: 70204

- Soft and comfortable earloop mask
- Ideal for sensitive skin
- Cool to wear
- Fluid resistant to 120 mm/Hg
- Complies to A\$4381 and ADA Guidelines

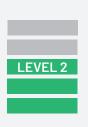



A34301		
Characteristics	Level 2 Barrier	Level 3 Barrier
Bacterial filtration efficiency (BFE), %	>98	>98
Differential pressure (ΔP), mm H²0/cm²	<5.0	<5.0
Resistance to penetration by synthetic blood, minimum pressure in mm Hg for pass result	120mm Hg	160mm Hg

AS4381

#### Applications

Level 2: For use in dentistry, where minimal blood droplet exposure may possibly occur.

Level 3: For all surgical procedures, major trauma first aid or in any area where the health care worker is at risk of blood or body fluid splash

#### Terms

Fluid Resistance: Resistance to penetration by synthetic blood under pressure showing no strike through on the inner facing. **BFE:** Bacterial Filtration is the measurement of the filtering ability of the mask when subjected to micron particles **Delta**  $\Delta P$ : The breathability of the mask. The lower the Delta  $\Delta P$ , the cooler and more comfortable the mask.

